
STANDARD Q

Dengue Duo

STANDARD™ Q Dengue Duo Test

PLEASE READ INSTRUCTIONS CAREFULLY BEFORE YOU PERFORM THE TEST

SD BIOSENSOR

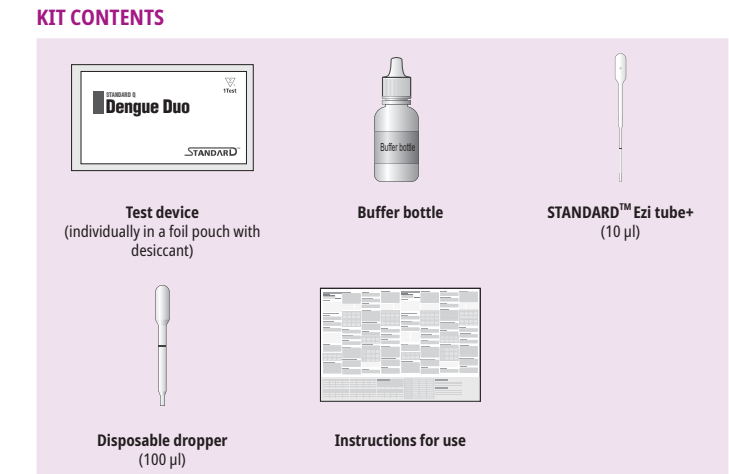MATERIALS REQUIRED BUT NOT PROVIDED

- Timer

- SPECIMEN COLLECTION**
- **Serum**
 - If serum in the plain tube is stored in a refrigerator at 2 ~ 8°C / 36 ~ 46°F, the specimen can be used for testing within 1 week after collection. Using the specimen in the long-term keeping more than 1 week can cause non-specific reaction. For prolonged storage, it should be at below -40°C / -40°F.
 - They should be brought to room temperature prior to use.
 - **Plasma**
 - If plasma in an anti-coagulant tube is stored in a refrigerator at 2 ~ 8°C / 36 ~ 46°F, the specimen can be used for testing within 1 week after collection. Using the specimen in the long-term keeping more than 1 week can cause non-specific reaction. For prolonged storage, it should be at below -40°C / -40°F.
 - They should be brought to room temperature prior to use.
 - **Whole blood**
 - **Capillary whole blood**
 - Capillary whole blood should be collected aseptically by fingrip.
 - Clean the area to be lanced with an alcohol swab.
 - Squeeze the end of the fingertip and pierce with a sterile lancet.
 - Collect the capillary whole blood to the black line of the STANDARD Ezi tube+ for the testing.
 - The capillary whole blood must be tested immediately after collection.
- [Venous Whole blood]**
- If venous whole blood in an anti-coagulant tube is stored in a refrigerator at 2 ~ 8°C / 36 ~ 46°F, the specimen can be used for testing within 1-2 days after collection.
 - Do not use hemolyzed blood specimens.

⚠

- Anticoagulants such as heparin or EDTA do not affect the test result.
- As known relevant interference, hemolytic specimen, rheumatoid factors-contained specimen and lipaemic, icteric specimen can lead to impair the test results.
- Use separate disposable materials for each specimen in order to avoid cross-contamination which can cause erroneous results.

PREPARATION AND TEST PROCEDURE

■ Preparation

Foil pouch

1

2

3

① Result window

② Specimen well

③ Buffer well

④ Yellow

⑤ Green

Test device

Desiccant

1. Carefully read the instruction for using the STANDARD Q Dengue Duo Test.

2. Look at the expiry date at the back of the foil pouch. Do not use the kit, if expiry date has passed.

3. Open the foil pouch, and check the test device and the Desiccant pack within the foil pouch.

- **Test Procedure**
[Dengue NS1 test device]
 - Using a sample collector (Disposable dropper (100 µl)), add 100 µl of the serum, plasma or whole blood to the specimen well of the test device.

100µl

Read After 15 mins

Do not read After 20 mins.

15-20 mins

[Dengue IgM/IgG test device]

1-1. Using a sample collector, add 10 µl of the serum, plasma or whole blood to the specimen well of the test device.

1-2. Add 3 drops (90 µl) of buffer to the buffer well of the test device.

PERFORMANCE CHARACTERISTICS				
■ The sensitivity and specificity of STANDARD Q Dengue Duo Test				
Total of 860 specimens were evaluated for sensitivity and specificity. The STANDARD Q Dengue Duo test kit got a high correlation with reference test (ELISA and RT-PCR).				
Reference	STANDARD Q Dengue Duo (NS1)			Total Result
	Positive	Negative		
RT-PCR	Positive	15	198	225
	Negative	222	225	
Total Result		186	423	
		Sensitivity: 183/198 (92.42%)		
		Specificity: 222/225 (98.7%)		

- Read the test results at 15-20 minutes. Do not read after 20 minutes.

⌚

Read After 15 mins

Do not read After 20 mins.

15-20 mins

⚠

Do not read test results after 20 minutes, it may give false results.

Test result		Example	
Dengue NS1	Negative		
	Positive		
	Invalid		
Dengue IgM/IgG	Negative		
	Dengue IgM Positive		
	Dengue IgG Positive		
	Dengue IgM/IgG Positive		
	Invalid		

- * The presence of any line no matter how faint it is, should be considered as a line formed.
- * Positive results should be considered in conjunction with the clinical history and other data available.
- * This test is for screening purposes. Confirmatory testing according to national guidelines is recommended to confirm the infection status.

EXPLANATION AND SUMMARY

- **Introduction**
Dengue viruses, transmitted by Aedes aegypti and Aedes albopictus mosquitoes, are widely distributed throughout the tropical and subtropical areas of the world. There are four known distinct serotypes of dengue virus (DEN-1, DEN-2, DEN-3 and DEN-4). Rapid and reliable tests for primary and secondary infections of Dengue are essential for patient management. An infected person experiences the acute symptoms of Dengue when there is a high level of the virus in the bloodstream. As the immune response fights the Dengue infection, the person's B cells begin producing IgMs and IgGs antibodies that are released in the blood and lymph fluid, where they recognize and neutralize the Dengue virus and viral molecules such as the Dengue non-structural protein 1 (NS1) antigen.

- **Intended use**
STANDARD Q Dengue Duo Test is an immunochromatographic assay for the detection of dengue NS1 antigens and IgM/IgG antibodies to dengue in human serum, plasma, or whole blood specimens. This test kit is for use in in-vitro diagnostics and for professional use.

- **Test principle**
STANDARD Q Dengue Duo Test Kit has two devices. For NS1 Ag device, the rapid test membrane is coated with an anti-Dengue NS1 on the test line. The specimen is added directly to the specimen well and interacts with monoclonal anti-Dengue NS1-gold on the conjugate pad. This specimen interacted with monoclonal anti-Dengue NS1-gold moves along membrane to the test line via capillary action to react with the anti-Dengue NS1. If NS1 antigen is present, a red line will appear at the test line. For IgM/IgG Ab device, monoclonal anti-human IgM and monoclonal anti-human IgG are immobilized at two individual test lines respectively (IgM/IgG line) on the membrane. Inactivated Dengue virus and monoclonal anti-Dengue-gold complex of pad release by adding buffer and react with anti-Dengue IgM or IgG. Human IgM and IgG in patient serum migrate and react with monoclonal anti-human IgM and IgG respectively on the membrane. If human anti-Dengue IgM or IgG exist in patient serum, the individual test line induced red colored band respectively, which means a positive test results. The red line at the control region should always appear if the assay is performed correctly.

KIT STORAGE AND STABILITY
Store the kit at 2 ~ 40°C / 36 ~ 104°F, out of direct sunlight. Kit materials are stable until the expiration date printed on the outer box. Do not freeze the kit.

WARNINGS AND PRECAUTIONS

- Do not re-use the test kit.
- Do not use the test kit if the pouch is damaged or the seal is broken.
- Do not use buffer bottle of another lot.
- Do not smoke, drink or eat while handling specimen.
- Wear personal protective equipment, such as gloves and lab coats when handling kit reagents. Wash hands thoroughly when afterwards.
- Clean up spills thoroughly using an appropriate disinfectant.
- Handle all specimens as if they contain infectious agents.
- Observe established precautions against microbiological hazards throughout testing procedures.
- Dispose of all specimens and materials used to perform the test as bio-hazard waste. Laboratory chemical and bio-hazard wastes must be handled and discarded in accordance with all local, state, and national regulations.
- Desiccant in foil pouch is to absorb moisture and keep humidity from affecting products. If the moisture indicating desiccant beads change from yellow to green, the test device in the pouch should be discarded.

LIMITATION OF TEST

- The test procedure, precautions and interpretation of results for this test must be followed strictly when testing.
- This test detects the presence of Dengue NS1 antigen and Dengue IgM/IgG antibodies in the specimen and should not be used as the sole criteria for the diagnosis of Dengue virus infection.
- Test results must be considered with other clinical data available to the physician.
- For more accuracy of immune status, additional follow-up testing using other laboratory methods is recommended.
- Neither the quantitative value nor the rate of Dengue NS1 Ag or anti-Dengue IgM/IgG concentration can be determined by this qualitative test.
- Failure to follow the test procedure and interpretation of test results may adversely affect test performance and/or produce invalid results.

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ANALYTICAL PERFORMANCE

- **Cross Reactivity**
For STANDARD Q Dengue Duo Test- Dengue NS1 Ag device, no cross-reactivity with Zika Virus, Chikungunya virus, Streptococcus pneumonia, Streptococcus pyogenes, Legionella pneumophila, Malaria plasmodium falciparum, Malaria plasmodium vivax, Influenza A (H1N1), Influenza A (H3N2), Influenza B, E.coli (K99), Salmonella (Typhi), Rota virus Type I, Rota virus Type II, Rota virus Type III, Rota virus Type IV, HBV antigen positive specimen, CMV Ab positive, vaccinated yellow fever virus, Japanese encephalitis virus, Hepatitis A virus, Rubella virus was observed.
For STANDARD Q Dengue Duo Test- Dengue IgM/IgG device, no cross-reactivity with HBV positive plasma, Chikungunya IgG positive plasma, Chikungunya IgM positive plasma, Zika IgM positive plasma, Syphilis Ab positive plasma, HIV-1 positive plasma, HIV-2 positive plasma, HBV antigen positive plasma, HCV positive plasma, CMV Ab positive specimen, vaccinated yellow fever virus, Japanese encephalitis, tick borne encephalitis antibody, Epstein Barr virus IgM, Parvovirus B19 antibody, and Varicella Zoster virus antibody was observed while it had cross reaction with Zika IgG, West Nile virus antibody, Chikungunya virus IgM and malaria Pf antibody.

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Manufactured by **SD Biosensor, Inc.**

Head office : C-4th&5th, 16, Deogyong-daero 1556beon-gil, Yeongtong-gu, Suwon-si, Gyeonggi-do, 16690, REPUBLIC OF KOREA

Manufacturing site : 74, Osongsangmyeong 4-ro, Osong-eup, Heungdeok-gu, Cheongju-si, Chungcheongbuk-do, 28161, REPUBLIC OF KOREA

Contact

e-mail: ts@sdbiosensor.com | phone: +82-80-970-9700 | website: www.sdbiosensor.com

SYMBOL

REF	Reference number	Caution	Use by	LOT	Batch code	Consult Instructions for Use	Do not re-use
IVD	In vitro Diagnostics	Note	Manufacturer		Date of manufacture	Contains Sufficient for <=> Tests	Keep away from sunlight
	Indicate that you should keep the product dry	To indicate the temperature limitations in which the transport package has to be kept and handled.		This product fulfills the requirements of the European Directive 98/79/EC.	Do not use if packaging is damaged		

BIBLIOGRAPHY

- Dengue guidelines for diagnosis, treatment, prevention and control. World Health Organization, New Edition 2009.
- Hills SC, Nimmanitong S, Nisalak A, Burke DS. Evidence that maternal dengue antibodies are important in the development of dengue hemorrhagic fever in infants. Am J Trop Med Hyg Jan, 38(2):411-419, 1988.
- Dengue haemorrhagic fever: Diagnosis, treatment, prevention and control, World Health Organization 2nd Edition, 1997.
- Ludolf D. et al., Serological differentiation of infections with dengue virus serotypes 1 to 4 by using recombinant antigens, J Clin Microbiol, 40(11):4317-4320, 2002.
- Matthew T. R. et al. Dengue virus pirates human platelets, Blood, 126(3):286-287, 2015.
- Guzman M. G. et al. Dengue: A continuing global threat, Nat Rev Microbiol, 8:57–516, 2010.

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MIT Promed Consulting GmbH